

# Start a case




1. Confirm surgeon profile
2. Start case from home screen
3. Pre-schedule patient case
4. Access settings

## Capturing images and video



1. Select the **video source**.
2. **Stream** video and images to a server
3. **Camera** to capture an image.
4. **Record/Stop** to start/stop recording live video
5. **Pause** recording
6. **Microphone** to mute/unmute audio recording
7. **Print Remaining Images** to print any images that have not been printed via auto-print
8. **Image/video count** to view captured images and videos in the Media Gallery


# Exporting and Printing media

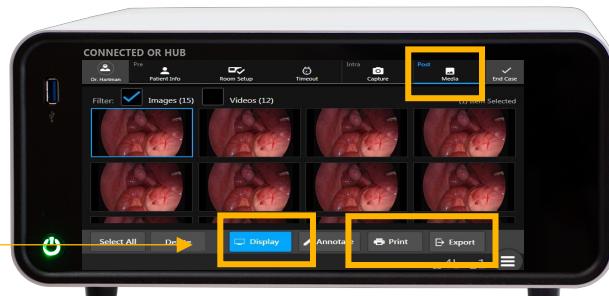

1. At the top of the Media Gallery, select **Images, Videos, or All**.
2. Select the images and videos to export or print
3. Press **Export or Print**
4. If **Exporting**, select the export destination (USB, Disc, DICOM, FTP, etc.)
5. If **Printing**, select print options: destination printer, image orientation, copies, and images per page.
6. Press **Export or Print**.
7. Press **x** to exit.

## Case icons

With DICOM or EMR Connectivity package, search information from PACS or EMR

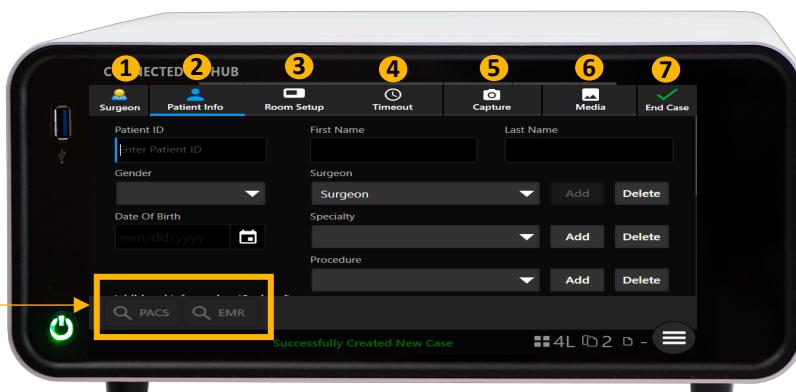

### 1. Surgeon

Select a surgeon profile

### 2. Patient Info

Enter patient information

### 3. Room Setup

Select the specialty, procedure type, and device preset settings

### 4. Timeout

Perform the surgical timeout checklist

### 5. Capture

Capture images and video

### 6. Media Gallery

View, annotate, and edit images

### 7. End Case

Finish the case for export or archival